


# 产品说明书

## CD11c Antibody

货号：ABB5206

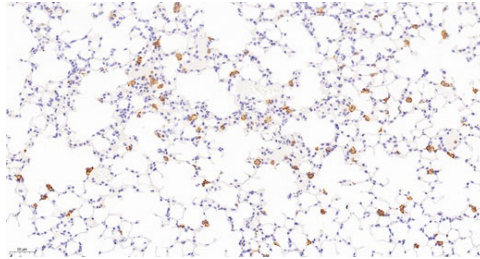

|           |                                                                                                                                                                                                                                        |
|-----------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 产品名称      | CD11c Antibody                                                                                                                                                                                                                         |
| 来源宿主      | Rabbit                                                                                                                                                                                                                                 |
| 反应种属      | Human                                                                                                                                                                                                                                  |
| 克隆类型      | Monoclonal                                                                                                                                                                                                                             |
| 同种型       | IgG                                                                                                                                                                                                                                    |
| 标记        | unconjugated                                                                                                                                                                                                                           |
| 纯化方式      | Affinity-chromatography                                                                                                                                                                                                                |
| 形式        | Liquid                                                                                                                                                                                                                                 |
| 存储溶液      | Rabbit IgG in phosphate buffered saline , pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.                                                                                                                                     |
| 存储方式      | Store at +4°C short term. Store at -20°C long term. Avoid freeze / thaw cycle.                                                                                                                                                         |
| 应用        | WB IHC-P/mlHC ICC/IF IP                                                                                                                                                                                                                |
| 使用方法      | WB 1:500~1:2000 IHC-P/mlHC 1:50~1:200 ICC/IF 1:50~1:200 IP 1:50                                                                                                                                                                        |
| 有效期       | one year                                                                                                                                                                                                                               |
| 别名        | p150/95 integrin alpha chain; Leu M5; SLEB6; ITGAX; CD11c; p150;Myeloid membrane antigen alpha subunit; Leukocyte surface antigen p150 95 alpha subunit; Integrin alpha-X;                                                             |
| 免疫原       | A synthesized peptide derived from human CD11c                                                                                                                                                                                         |
| SwissProt | P20702                                                                                                                                                                                                                                 |
| 分子量       | 128kDa                                                                                                                                                                                                                                 |
| 产品介绍      | Integrin alpha-X/beta-2 is a receptor for fibrinogen. It recognizes the sequence G-P-R in fibrinogen. It mediates cell-cell interaction during inflammatory responses. It is especially important in monocyte adhesion and chemotaxis. |

For research use only !

**慧蓝生物 竭诚为您服务**

地址：上海市浦东新区周浦镇天雄路588弄，电话：19101712317， 邮箱：469340997@qq.com