


产品说明书

CD31 Antibody

货号: ABB5338

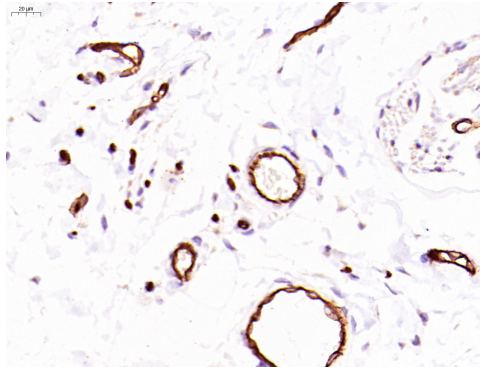

产品名称	CD31 Antibody
来源宿主	Rabbit
反应种属	Human
克隆类型	Monoclonal
同种型	IgG
标记	unconjugated
纯化方式	Affinity-chromatography
形式	Liquid
存储溶液	Rabbit IgG in phosphate buffered saline , pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.
存储方式	Store at +4°C short term. Store at -20°C long term. Avoid freeze / thaw cycle.
应用	WB IHC-P/mlHC IP FC
使用方法	WB 1:500~1:2000 IHC-P/mlHC 1:50~1:200 IP 1:50 FC 1:50
有效期	one year
别名	adhesion molecule; CD31 antigen; CD31/EndoCAM; EndoCAM; PECA1; PECAM-1; PECAM-1, CD31/EndoCAM; PECAM1; Platelet endothelial cell adhesion molecule;
免疫原	A synthesized peptide derived from human CD31
SwissProt	P16284
分子量	125kDa
产品介绍	PECAM-1 a cell adhesion molecule expressed on platelets and at endothelial cell intercellular junctions. Plays a role in modulation of integrin-mediated cell adhesion, transendothelial migration, angiogenesis, apoptosis, cell migration, negative regulatio

For research use only !

慧蓝生物 竭诚为您服务

地址: 上海市浦东新区周浦镇天雄路588弄, 电话: 19101712317, 邮箱: 469340997@qq.com