


# 产品说明书

ly6g

货号: ABB00006

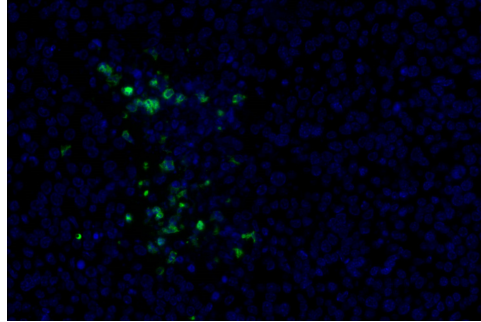

产品名称	ly6g
来源宿主	Rabbit
反应种属	Mouse
克隆类型	Rabbit monoclonal
克隆号	R9-135
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified
形式	Liquid
存储溶液	PBS (pH7.4),0.05% BSA,40% Glycerol. Preservative: 0.05% Sodium Azide.Stable for 12 months.
存储方式	Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles.Stable for 12 months.
应用	Flow Cyt, IHC-P,mIHC
使用方法	Flow Cyt 1:50; IHC-P 1:50-1:200; mIHC 1:50-1:200
有效期	one year
别名	GR-1,GR1,Lymphocyte antigen 6G
免疫原	Recombinant fragment.
SwissProt	P35461
细胞定位	Cell membrane
产品介绍	Ly6G is a GPI-anchored protein, that is also known as the myeloid differentiation antigen Gr1. The antigen is transiently expressed on monocytes in the bone marrow. The level of antigen expression in the bone marrow directly correlates with granulocyte differentiation and maturation. Ly6G is expressed predominantly on neutrophils, also in a subset of eosinophils, differentiating pre-monocytes, and plasmacytoid dendritic cells.

For research use only!

**慧蓝生物 竭诚为您服务**

地址: 上海市浦东新区周浦镇天雄路588弄, 电话: 19101712317, 邮箱: 469340997@qq.com