


# 产品说明书

## Vimentin Antibody

货号: ABB5134

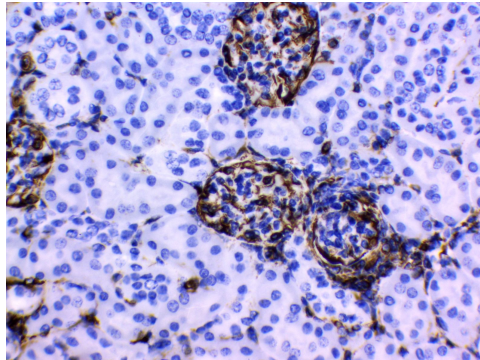

产品名称	Vimentin Antibody
来源宿主	Rabbit
反应种属	Human Mouse Rat
克隆类型	Monoclonal
同种型	IgG
标记	unconjugated
纯化方式	Affinity-chromatography
形式	Liquid
存储溶液	Rabbit IgG in phosphate buffered saline , pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.
存储方式	Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles.
应用	WB IHC-P/mlHC ICC/IF FC
使用方法	WB 1:500~1:2000 IHC-P/mlHC 1:50~1:200 ICC/IF 1:50~1:200 FC 1:30
有效期	one year
别名	FLJ36605; VIM; VIME; Vimentin
免疫原	A synthesized peptide derived from human Vimentin
SwissProt	P08670
分子量	54kDa
产品介绍	Vimentins are class-III intermediate filaments found in various non-epithelial cells, especially mesenchymal cells. Vimentin is attached to the nucleus, endoplasmic reticulum, and mitochondria, either laterally or terminally. Involved with LARP6 in the stabilization of type I collagen mRNAs for CO1A1 and CO1A2.

For research use only !

**慧蓝生物 竭诚为您服务**

地址: 上海市浦东新区周浦镇天雄路588弄, 电话: 19101712317, 邮箱: 469340997@qq.com