


## 产品说明书

## FOXP3

货号: ABB00009

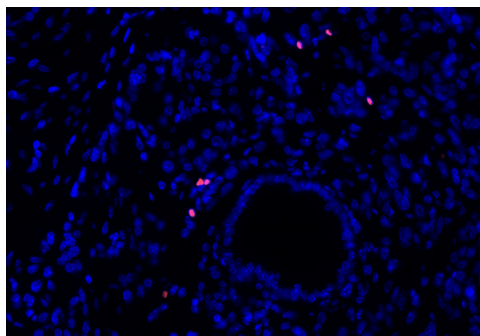

产品名称	FOXP3
来源宿主	Rabbit
反应种属	Human, Mouse, Rat
克隆类型	Rabbit monoclonal
克隆号	07-47
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified
形式	Liquid
存储溶液	PBS (pH7.4), 0.05% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide. Stable for 12 months.
存储方式	Store at -20°C. Avoid repeated freeze / thaw cycles. Stable for 12 months.
应用	WB, IHC-P, mIHC
使用方法	WB: 1:1000; IHC-P: 1:50-1:200, mIHC 1:50-1:200
有效期	one year
别名	Forkhead box P3; Forkhead box protein P3; foxp3
免疫原	Recombinant protein within human FOXP3 aa 101-370
SwissProt	Q9BZS1
细胞定位	Nucleus, Cytoplasm.
分子量	47kDa
产品介绍	FOXP3 (forkhead box P3), also known as scurfin, is a protein involved in immune system responses. A member of the FOX protein family, FOXP3 appears to function as a master regulator of the regulatory pathway in the development and function of regulatory T cells. Regulatory T cells generally turn the immune response down. In cancer, an excess of regulatory T cell activity can prevent the immune system from destroying cancer cells. In autoimmune disease, a deficiency of regulatory T cell activity can allow other autoimmune cells to attack the body's own tissues.

For research use only!

慧蓝生物 竭诚为您服务

地址: 上海市浦东新区周浦镇天雄路588弄, 电话: 19101712317, 邮箱: 469340997@qq.com