


产品说明书

CD68

货号：ABB00011

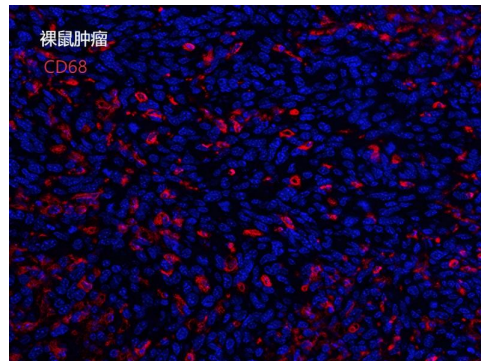

产品名称	CD68
来源宿主	Rabbit
反应种属	Mouse,Rat
克隆类型	Rabbit monoclonal
克隆号	22-85
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified
形式	Liquid
存储溶液	PBS (pH7.4),0.05% BSA,40% Glycerol. Preservative: 0.05% Sodium Azide.Stable for 12 months.
存储方式	Store at +4℃ after thawing. Aliquot store at -20℃. Avoid repeated freeze / thaw cycles.Stable for 12 months.
应用	WB,IF-Cell,IHC-P,mIHC IHC-Fr,IF-Tissue
使用方法	WB:1:1000;IF-Cell: 1:100;IHC-P:1:100-1:500;mIHC:1:100-1:500; IHC-Fr:1:50-1:200;IF-Tissue:1: 50-1: 200
有效期	one year
别名	CD68 antigen; CD68 molecule; Gp110; LAMP4; Macrophage antigen CD68 (microsialin); Macrosialin; SCARD1; Scavenger receptor class D member 1;
免疫原	Recombinant protein within mouse CD68.
SwissProt	P31996
细胞定位	Cell membrane and Endosome membrane. Lysosome membrane.
分子量	37kDa
产品介绍	Could play a role in phagocytic activities of tissue macrophages, both in intracellular lysosomal metabolism and extracellular cell-cell and cell-pathogen interactions. Binds to tissue- and organ-specific lectins or selectins, allowing homing of macrophage subsets to particular sites. Rapid recirculation of CD68 from endosomes and lysosomes to the plasma membrane may allow macrophages to crawl over selectin-bearing substrates or other cells.

For research use only !

慧蓝生物 竭诚为您服务

地址：上海市浦东新区周浦镇天雄路588弄，电话：19101712317，邮箱：469340997@qq.com