


# 产品说明书

CD8α

货号：ABB00014

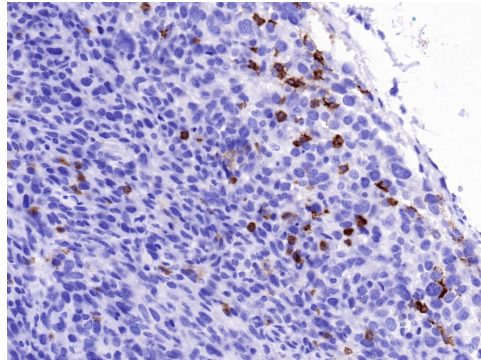

产品名称	CD8α
来源宿主	Rabbit
反应种属	Mouse
克隆类型	Rabbit monoclonal
克隆号	42-52
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified
形式	Liquid
存储溶液	PBS (pH7.4), 0.05% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide. Stable for 12 months.
存储方式	Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles. Stable for 12 months.
应用	IHC-P, mIHC, WB, FC(Extra), IF
使用方法	IHC-P: 1: 100-1: 500; mIHC: 1: 100-1: 500; WB: 1: 1000; FC(Extra): 1: 50; IF: 1: 100-1: 200
有效期	one year
别名	CD8a, Lyt-2, T-cell surface glycoprotein Lyt-2
免疫原	Synthetic Peptide.
SwissProt	P01731
细胞定位	Cell membrane
分子量	27kDa
产品介绍	CD8 is a transmembrane glycoprotein that is predominantly expressed on the surface of cytotoxic T cells, and can also be found on natural killer cells, cortical thymocytes, and dendritic cells. CD8 serves as a co-receptor for the T cell receptor (TCR). Both CD8 and TCR recognize antigens displayed by an antigen presenting cell (APC) in the context of class I MHC molecules. CD8 plays a role in T cell development and activation of mature T cells. This antibody is raised against 23-183 aa of human CD8a, the alpha chain of CD8 molecule.

For research use only !

**慧蓝生物 竭诚为您服务**

地址：上海市浦东新区周浦镇天雄路588弄，电话：19101712317，邮箱：469340997@qq.com