


产品说明书

CD11b

货号：ABB00015

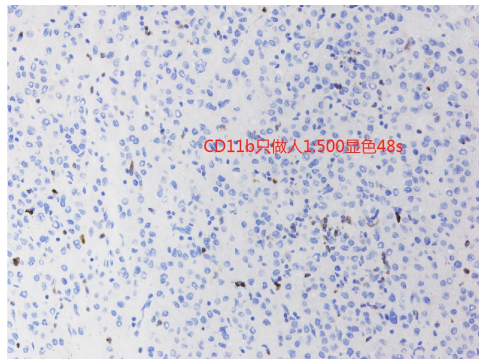

产品名称	CD11b
来源宿主	Rabbit
反应种属	Human
克隆类型	Rabbit monoclonal
克隆号	58-44
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified
形式	Liquid
存储溶液	PBS (pH7.4), 0.05% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide. Stable for 12 months.
存储方式	Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles. Stable for 12 months.
应用	IHC-P, mIHC, ICC, WB, IP, FC(Intra)
使用方法	IHC-P: 1:100-1:500; mIHC: 1:100-1:500; ICC: 1:100-1: 500; WB: 1:200-1: 1000; IP: 1: 25
有效期	one year
别名	CD11B, CD11B/Integrin Alpha M, CR 3 alpha chain, CR3A, Integrin alpha M, ITGAM, MAC 1, MAC1A, MO1A, Neutrophil adherence receptor, SLEB6
免疫原	Synthetic Peptide
SwissProt	P11215
细胞定位	Cell membrane
分子量	127kDa
产品介绍	CD11b, also known as Integrated alpha-m, a transgender protein, can form an heterododorous composed of α and β subunit. It is a common bone marrow mark (neutrophil, monocyte, macrophage, and small gel cells) and NK (natural kill cells) antigens. It can be used to distinguish between acute granulocyte deficiency (CD11B+, CD117-) and acute early early early elastic cell leukemia (CD11B-, CD117+).

For research use only !

慧蓝生物 竭诚为您服务

地址：上海市浦东新区周浦镇天雄路588弄，电话：19101712317，邮箱：469340997@qq.com