


产品说明书

Mannose Receptor/CD206 Antibody

货号：ABB00004

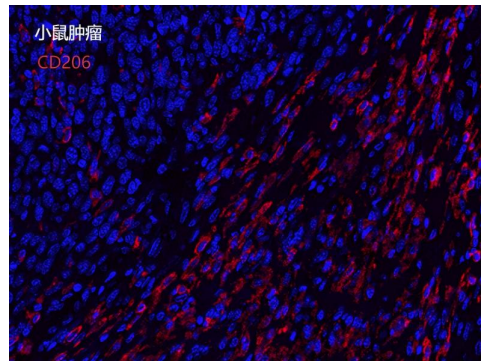

产品名称	Mannose Receptor/CD206 Antibody
来源宿主	Rabbit
反应种属	Mouse,Rat
克隆类型	Rabbit monoclonal
克隆号	07-76
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified
形式	Liquid
存储溶液	PBS (pH7.4),0.05% BSA,40% Glycerol. Preservative: 0.05% Sodium Azide.Stable for 12 months.
存储方式	Store at +4℃ after thawing. Aliquot store at -20℃. Avoid repeated freeze / thaw cycles.Stable for 12 months.
应用	IHC-P, mlHC, IF
使用方法	IHC-P: 1: 50-1: 200; mlHC: 1: 50-1: 200; IF: 1: 50-1: 100
有效期	one year
别名	CD206, Macrophage mannose receptor 1 , MRC1, Mannose receptor C type 1
免疫原	Recombinant protein within mouse Mannose Receptor(CD206) aa 1-1,409.
SwissProt	Q61830
细胞定位	Endosome membrane, Cell membrane
分子量	165kDa
产品介绍	Mediates the endocytosis of glycoproteins by macrophages. Binds both sulfated and non-sulfated polysaccharide chains. Acts as phagocytic receptor for bacteria, fungi and other pathogens.

For research use only !

慧蓝生物 竭诚为您服务

地址：上海市浦东新区周浦镇天雄路588弄，电话：19101712317，邮箱：469340997@qq.com