


产品说明书

CD86

货号：ABB00019

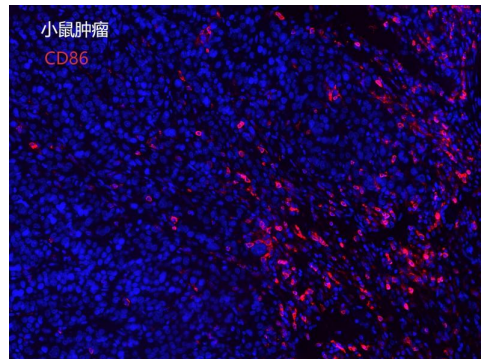

产品名称	CD86
来源宿主	Rabbit
反应种属	Mouse
克隆类型	Rabbit monoclonal
克隆号	00-48
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified
形式	Liquid
存储溶液	PBS (pH7.4), 0.05% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide. Stable for 12 months.
存储方式	Store at -20°C. Avoid repeated freeze / thaw cycles. Stable for 12 months.
应用	WB, IHC-P, mIHC, IF, IP, ELISA
使用方法	IHC-P 1:200, mIHC 1:200, WB 1:500-1:5000, IF 1:200, ELISA 1:5000-1:20000, IP 1:50-200
有效期	one year
别名	B70; B7-2; B7.2; LAB72; CD28LG2; T-lymphocyte activation antigen CD86
免疫原	Synthetic peptide.
Gene ID	1
SwissProt	P42082
细胞定位	Membrane.
分子量	60-80 kDa
产品介绍	Receptor involved in the costimulatory signal essential for T-lymphocyte proliferation and interleukin-2 production, by binding CD28 or CTLA-4. May play a critical role in the early events of T-cell activation and costimulation of naive T-cells, such as deciding between immunity and anergy that is made by T-cells within 24 hours after activation. Isoform 2 interferes with the formation of CD86 clusters, and thus acts as a negative regulator of T-cell activation.

For research use only !

慧蓝生物 竭诚为您服务

地址：上海市浦东新区周浦镇天雄路588弄，电话：19101712317，邮箱：469340997@qq.com