


产品说明书

CD31 (Y0350) Rabbit mAb

货号: HL18207

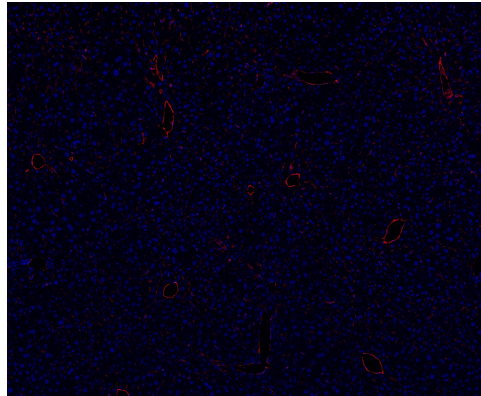

产品名称	CD31 (Y0350) Rabbit mAb
来源宿主	Rabbit
反应种属	Human, Mouse, Rat,
克隆类型	Rabbit monoclonal
克隆号	Y0350
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified.
形式	Liquid
存储溶液	PBS (pH7.4), 50% Glycerol, 0.05% BSA, 0.05% Proclin 300
存储方式	Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles.
应用	WB, IHC, mlHC, IF
使用方法	IHC 1:500-1:1000; mlHC 1:500-1:1000; WB 1:500-1:2500; IF 1:100-1:500
有效期	one year
别名	Platelet endothelial cell adhesion molecule; PECAM-1; EndoCAM; GPIIA; PECA1; CD antigen CD31;
SwissProt	P16284
细胞定位	Membrane
分子量	130kD
产品介绍	platelet and endothelial cell adhesion molecule 1(PECAM1) Homo sapiens The protein encoded by this gene is found on the surface of platelets, monocytes, neutrophils, and some types of T-cells, and makes up a large portion of endothelial cell intercellular junctions. The encoded protein is a member of the immunoglobulin superfamily and is likely involved in leukocyte migration, angiogenesis, and integrin activation. [provided by RefSeq, May 2010],

For research use only !

慧蓝生物 竭诚为您服务

地址: 上海市浦东新区周浦镇天雄路588弄, 电话: 19101712317, 邮箱: 469340997@qq.com