


# 产品说明书

GFAP (Y0298) Rabbit mAb

货号: HL18172

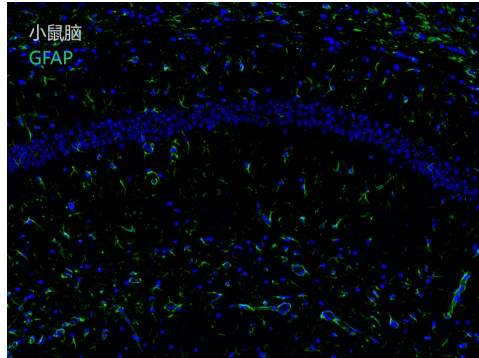

|           |   |
|-----------|---|
| 产品名称      | GFAP (Y0298) Rabbit mAb   |
| 来源宿主      | Rabbit  |
| 反应种属      | Human, Mouse, Rat,  |
| 克隆类型      | Rabbit monoclonal   |
| 克隆号       | Y0298   |
| 同种型       | IgG   |
| 标记        | unconjugated  |
| 纯化方式      | Protein A affinity purified.  |
| 形式        | Liquid  |
| 存储溶液      | PBS (pH7.4), 50% Glycerol, 0.05% BSA, 0.05% Proclin 300   |
| 存储方式      | Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles.   |
| 应用        | WB, IHC, mlHC, IF   |
| 使用方法      | IHC 1:100-1:500; mlHC 1:100-1:500; WB 1:500-1:2500; IF 1:50-1:200   |
| 有效期       | one year  |
| 别名        | GFAP; Glial fibrillary acidic protein; GFAP   |
| SwissProt | P14136  |
| 细胞定位      | Cytoplasm   |
| 分子量       | 50kD  |
| 产品介绍      | This gene encodes one of the major intermediate filament proteins of mature astrocytes. It is used as a marker to distinguish astrocytes from other glial cells during development. Mutations in this gene cause Alexander disease, a rare disorder of astrocytes in the central nervous system. Alternative splicing results in multiple transcript variants encoding distinct isoforms. [provided by RefSeq, Oct 2008], |

For research use only!

慧蓝生物 竭诚为您服务

地址: 上海市浦东新区周浦镇天雄路588弄, 电话: 19101712317, 邮箱: 469340997@qq.com