


产品说明书

α -SMA Mouse mAb

货号: HL06566

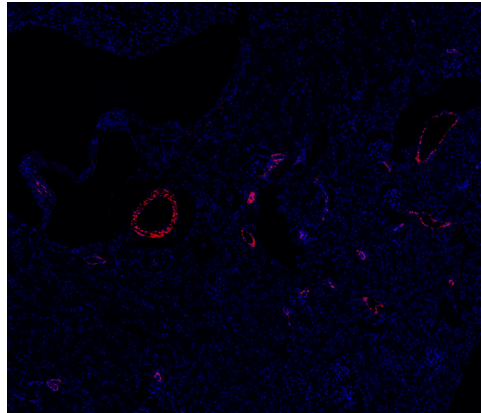

产品名称	α -SMA Mouse mAb
来源宿主	Mouse
反应种属	Human Mouse Rat
克隆类型	Mouse monoclonal
同种型	IgG
标记	unconjugated
纯化方式	Protein A affinity purified.
形式	Liquid
存储溶液	PBS (pH7.4), 0.05% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide.
存储方式	Store at +4°C after thawing. Aliquot store at -20°C. Avoid repeated freeze / thaw cycles.
应用	IHC-P, mlHC
使用方法	IHC-P 1:100-1:1000; mlHC 1:100-1:1000
有效期	one year
别名	ACTA2/ACTG2, Actin, smooth muscle
免疫原	Synthesized peptide derived from human Actin, smooth muscle (SMA) AA range: 1-100
SwissProt	P62736
细胞定位	Cytoplasmic
产品介绍	The protein encoded by this gene belongs to the actin family of proteins, which are highly conserved proteins that play a role in cell motility, structure and integrity. Alpha, beta and gamma actin isoforms have been identified, with alpha actins being a major constituent of the contractile apparatus, while beta and gamma actins are involved in the regulation of cell motility. This actin is an alpha actin that is found in skeletal muscle. Defects in this gene cause aortic aneurysm familial thoracic type 6. Multiple alternatively spliced variants, encoding the same protein, have been identified.

For research use only !

慧蓝生物 竭诚为您服务

地址: 上海市浦东新区周浦镇天雄路588弄, 电话: 19101712317, 邮箱: 469340997@qq.com